

5 MEGAPIXEL HR

AM7515MT8A Dino-Lite Edge - Polarisasi



Dipelopori dengan built-in lampu koaksial, AM7515MT8A yang menjawab kebutuhan pengamatan brightfield dengan kualitas gambar yang superior pada pembesaran 700 ~ 900x. Melalui penggunaan Fleksibel LED Control (FLC), AM7515MT8A tidak hanya dapat beralih antara brightfield dan darkfield tetapi juga dapat menggabungkan secara bersama untuk lebih banyak kemungkinan pada menyoroti objek secara detail. AM7515MT8A menyediakan Automatic Magnification Readings (AMR) yang dirancang untuk mengoptimalkan pengukuran dengan efisiensi dan akurasi.


The AM7515MT8A menyediakan cahaya koaksial (brightfield iluminasi) yang mirip dengan fungsi dengan menghubungkan N3C-A.

Cahaya melewati lensa untuk membentuk cahaya aksial yang ideal untuk melihat kasar dan tidak rata pada permukaan, seperti goresan, penyok, karakter yang dicetak, permukaan cermin, dan benda sangat reflektif.


Ini ideal untuk meningkatkan resolusi gambar di bawah pembesaran tinggi.

IPHONE PANEL

Coaxial Light (Bright field illumination)

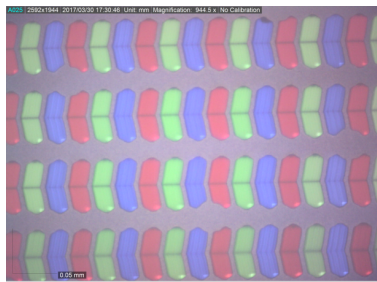

Dark Illumination

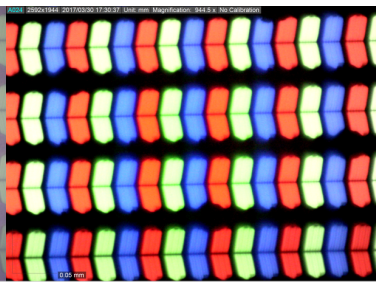

Pertanyaan **contact support (Klik)**

Warranty : 1 Tahun, Made in : Taiwan

Optional :

Basic Stand

Profesional Stand



SPECIFICATIONS

Special feature	Auto Magnification Reading, FLC (Switch LED), Brightfield / Darkfield
ESD Safe	Yes
Measurement functionality	Yes
Housing	Aluminum
Cap Changeable	Yes (AD / EDGE series)
Polarizer	Yes
Number of LEDs	8 (Switch FLC)
Light / LED type	White, Brightfield / Darkfield
Working distance	Long
Magnification	700~900x
Resolution	5 MP (2592 x 1944)
Connection	USB 2.0

ADDITIONAL INFO

Model	AM7515MT8A Dino-Lite Edge
Interface	USB 2.0
Product Resolution	5M pixels (2592x1944)
Magnification	700x~900x
Frame Rate	10fps in 5MP/3MP/2MP, MJPEG 25fps in 1.3MP, MJPEG 30fps in VGA, MJPEG
Lighting	8 white LEDs
Microtouch	Yes
Operating System Supported	Windows 10, 8, 7, Vista, XP (Windows OS Only)
Calibration Function	Yes
Measurement Function	Yes
Unit Dimension	10.7cm (H) x 3.2cm (D)
Unit Weight	139g

Information about working distance and field of view

M	WD	FOV (x)	FOV (y)	DOF
700	6.6	0.6	0.4	0.09
800	6.4	0.5	0.4	0.08
900	6.3	0.4	0.3	0.07

M = magnification rate

WD = working distance (without front cap)

FOV = field of view

DOF= depth of field

Unit = mm